


Imagine... consistently superior  
cytology preparation



**Tissue-Tek® Cyto-Tek® 2500**

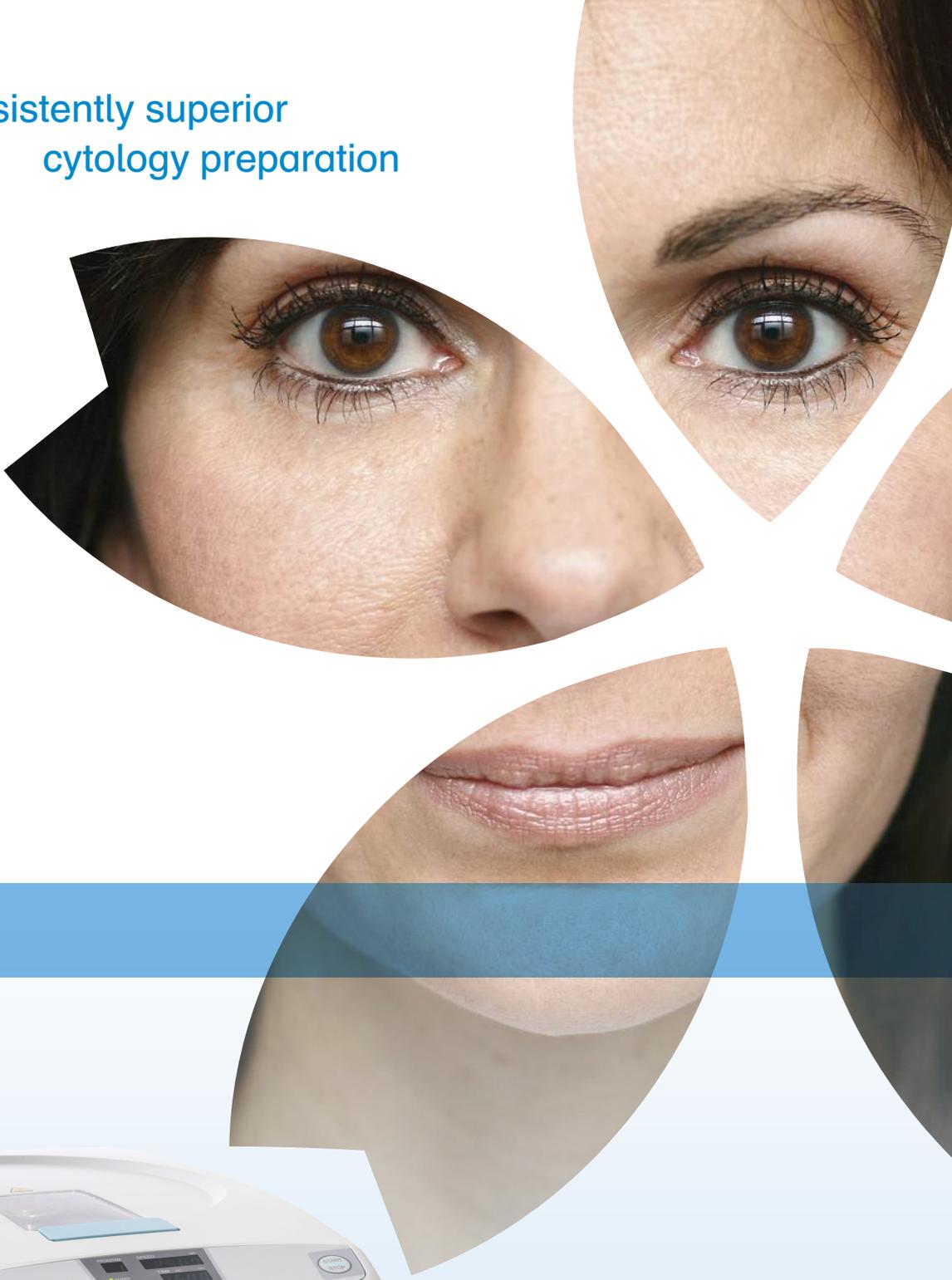

Cytocentrifuge

# TISSUE-TEK® CYTO-TEK® 2500 CYTOCENTRIFUGE

The Cyto-Tek® 2500 Centrifuge provides optimal cell recovery and advanced preservation of cellular structure in non-gynecological monolayer slide preparation.

The 2500 was designed with enhanced features focusing on user safety, including the introduction of a new sealable rotor housing with lid which significantly reduces aerosol exposure.

The housing is completely removable and transportable, providing added convenience. Technicians now have the ability to load and unload specimens in vented areas away from the centrifuge, under a fume hood.


The Cyto-Tek® 2500 Centrifuge has a sleek and modern design and provides the superior clinical performance users have come to expect from Sakura's Cyto-Tek® Centrifuges.

## Dimensions

W x D x H 375 x 490 x 225 mm  
Weight 11.7 kg

## Features

Capacity 12  
RPM 200 ~ 2.500 by 10rpm increments  
Program Memory 30  
Autoclavable Yes (Guard Bowl)

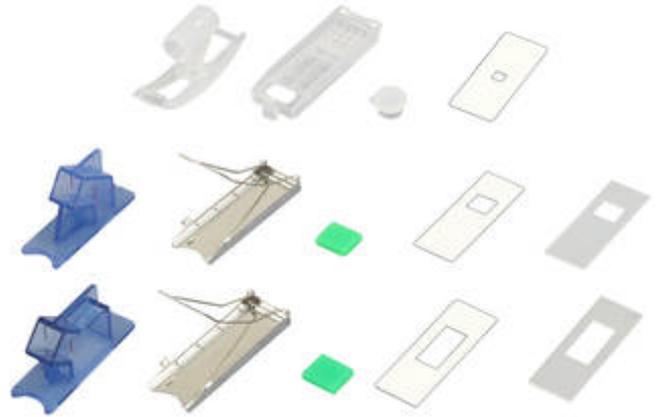

## ITEMS & OPTIONS



- |             |  |                   |                                       |
|-------------|--|-------------------|---------------------------------------|
| <b>4300</b> | Cyto-Tek®2500 Cytocentrifuge, 100-240V; 50/60 Hz | <b>F60-900-00</b> | Secondary box for Cyto-Tek®2500, 230V |
| <b>4301</b> | Chamber Kit 1mL; 50/pcs                          | <b>4338</b>       | Cyto-Tek®6 mL Filter Paper; 25/pcs    |
| <b>4326</b> | Cyto-Tek®6/12 mL Base Holder; 12/pcs             | <b>4339</b>       | Cyto-Tek®12 mL Filter Paper; 100/pcs  |
| <b>4327</b> | Cyto-Tek®6 mL Gasket; 25/pcs                     | <b>8704</b>       | Giemsa; 2.5 L                         |
| <b>4328</b> | Cyto-Tek®12 mL Fluid Chamber; 12/pcs             | <b>8705</b>       | Giemsa; 5 L                           |
| <b>4329</b> | Cyto-Tek®1 mL Fluid Chamber; 200/pcs             | <b>8716</b>       | Papanicolaou 2A OG6; 2.5 L            |
| <b>4331</b> | Cyto-Tek®6 mL Fluid Chamber; 12/pcs              | <b>8717</b>       | Papanicolaou 2A OG6; 5 L              |
| <b>4333</b> | Cyto-Tek®1 mL Filter Paper; 200/pcs              | <b>8719</b>       | Papanicolaou 2B Orange II; 2.5 L      |
| <b>4334</b> | Cyto-Tek®6/12 mL Chamber Cap; 25/pcs             | <b>8720</b>       | Papanicolaou 2B Orange II; 5 L        |
| <b>4336</b> | Cyto-Tek®1 mL Base Holder; 200/pcs               | <b>8722</b>       | Papanicolaou 3B EA 50; 2.5 L          |
| <b>4337</b> | Cyto-Tek®12 mL Gasket; 25/pcs                    | <b>8723</b>       | Papanicolaou 3B EA 50; 5 L            |
|             |  | <b>8731</b>       | Cervix Spray Fixative; 125 mL         |

**First we understand.  
Then we innovate.**



Sakura Finetek Europe B.V.  
The Netherlands  
Tel: +31 (0)88 592 00 00  
Fax: +31 (0)88 592 00 01  
Sakura@sakura.eu  
www.sakura.eu