

Slide Ultraviolet Printer G363UV




Size: (L) 457* (W) 275* (H) 517mm

Slide Ultraviolet Printer G363UV, using the latest UV laser printing technology, can clearly print high-resolution patterns and information on the slide coating.

Print speed:5s/piece

Dot per inch:2500dpi

Load/Collection capacity:up to 150 pieces of slides.

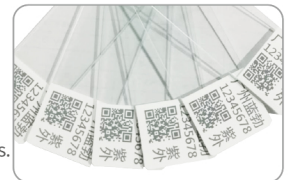

/Technical Parameters/

- **Latest UV laser printing technology**

Small laser spot and higher fineness can clearly print out the extremely high-resolution QR code and related information.

- **Zero cost**

No ribbon or ink jet. Slide without adding laser powder (radium carving powder). Except slide, no other consumable cost.

- **HD LCD screen**

10 inch HD LCD touch screen.

Language: English and others as customer request.

- **Built-in operating system without external connection**

Stand-alone work with built-in computer .

- **Clear print, long-term preservation.**

Print Texts, Symbol, Digits, Graphics and QR code.

Printed contents are corrosion resistant exposed to xylene and alcohol.

- **Multiple information input methods**

Input through slide ultraviolet printer, Tissue Cassette Scanner, remote network, LIS/HIS system etc.

- **Light absorbing material**

Using light absorbing material to avoid light reflection and light pollution effectively.

